



04450A

Please read the following instructions and warnings carefully.

Keep this instruction manual for future reference.

No Tools Required

Requires 1 "AA" Battery (not included)

**For Technical Support Call 1-800-268-6237 or
E-Mail at customerservice@summerinfant.com**

For use with children from birth and up.

 **WARNING:**
Adult assembly required.
Keep small parts away from children when assembling.

INSTRUCTIONS

GENTLE VIBRATIONS MASSAGER & TOOTHBRUSH

Thank you for purchasing the Gentle Vibrations Massager & Toothbrush by Summer Infant. The toothbrush includes a battery-operated handle, 2 adaptable toothbrush heads, and toothbrush stand with suction cup mounting. The soothing vibrating heads gently brush baby's gums and teeth. The first adaptable head has a mirror for viewing baby's new teeth on one side and soft toothbrush bristles for toddler on the other. The second adaptable head has a gum massager on one side and a rubber infant toothbrush on the other side. The suction cup on base allows for upright storage.


BATTERY WARNING:

- Do not mix old and new batteries.
- Do not mix alkaline, standard, or rechargeable batteries.
- Remove batteries when product is stored for long period of time or when batteries are exhausted.
- Dispose of exhausted batteries properly.

CAUTION

This product is not a toy. Adult supervision required.
Store away from child's reach.

FEATURES & COMPONENTS:

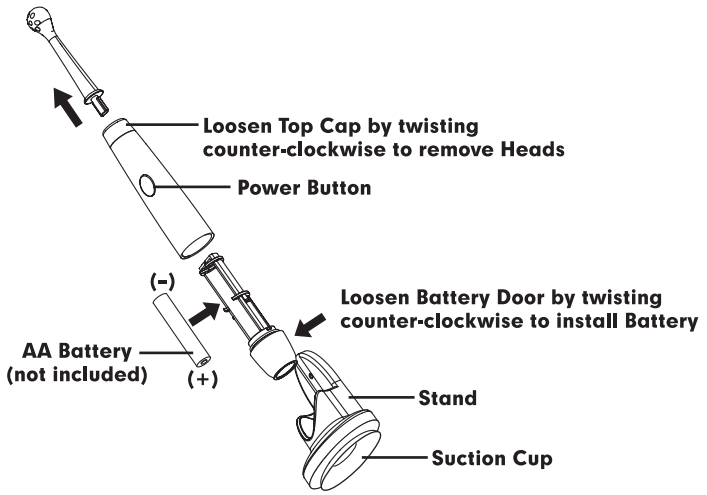

BATTERY INSTALLATION

This product requires 1 "AA" battery (not included). Twist battery cap counter-clockwise at base of handle and insert battery. Replace battery cap and twist battery cap clockwise to tighten.

ADAPTABLE HEADS

This unit comes with two adaptable heads:

- Head with mirror for viewing baby's new teeth on one side and soft toothbrush bristles for toddlers on the other (pre-installed)
- Head with gum massager on one side and infant toothbrush on the other (in stand)

To change the heads, loosen top cap by twisting top cap counter-clockwise and remove installed adaptable head by pulling head up. Plug alternate head into top cap and twist cap clockwise until notch on top cap lines up with LOCK notch on main body.

Additional toothbrush heads may be purchased separately. Feel free to contact us at the below phone number or email if you should need to order replacement or additional products.

We at Summer Infant stand behind all of our products. If you are not completely satisfied or have any questions, please contact our Customer Service Department at 1-800-268-6237 or e-mail us at customerservice@summerinfant.com.